

Anti-CD147 Monoclonal Antibody

Purified antibody

(BSG basigin (Ok blood group))

code Number : CD147-Mo-MY031 (RRID : AB_2753114)

Clone name : FRT35

Other name : OK; 5F7; TCSF; EMPRIN; EMMPRIN; SLC7A11

Size : 0.1mg (affinity-purified, proteinG)

Formulation : liquid. (1.1 mg / ml in PBS with 0.03% ProClin)

Storage : Store at cool temp (2-10°C).

The antibody can be stored at 2-10°C. This antibody can also be aliquotted and stored at -80°C for long-term storage. Avoid repeated freeze-thaw cycles. Non-hazardous. No MSDS required.


Species : mouse, monoclonal

Isotype : IgG1, κ


Antigen : NCI-H322/Human

Specificity : human (others not tested)

This antibody is specific to CD147.

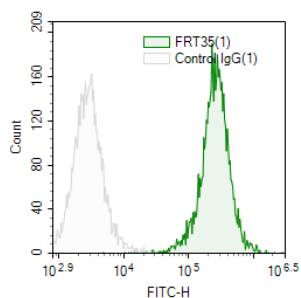

cell line:KG1A
2nd antibody:anti-mouse IgG FITC

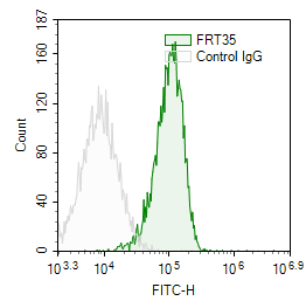

cell line:IM9
2nd antibody:anti-mouse IgG FITC

ご注意 : 本商品には 0.03 % の ProClin が入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD

Applications : Flow Cytometry (FC), ImmunoPrecipitation (IP)

Research Use : For research use only, not for use in diagnostic procedures.

Remark : FC recommended concentration

(cell conc. 5×10^5 cell / test & antibody conc. 0.5 μg / test)

Reference :

ご注意：本商品には 0.03 % の ProClin が入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD