

## *Anti-CD81 Monoclonal Antibody*

### *Purified antibody*

*(TAPA1, Tspan-28)*

#### Product Information:

---

1) Catalog #:	MSFR100760
Item Name:	CD81 mAb FRT6 (purified) 100ug
Size :	100µg (affinity-purified, protein G)
Species :	Mouse
Product Code :	CD81-Mo-MY005
RRID :	AB_2722744
2) Catalog #:	MSFR107350
Item Name:	CD81 mAb FRT6 (purified) 1 mg
Size :	1mg (affinity-purified, protein G)
Species :	Mouse
Product Code :	CD81-Mo-MY005
RRID :	AB_2722744

---

**Clone name :** FRT6

**Other name :** TAPA1, Tspan-28

**Formulation :** liquid. ( 1 mg / ml in PBS )

**Storage :** Store at 4°C. This antibody can also be aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thaw cycles. Non-hazardous. No MSDS required.


**Species :** mouse, monoclonal

**Isotype :** IgG2a, κ


**Antigen :** SKOV3/Human

**Specificity :** human (others not tested)

This antibody is specific to CD81.

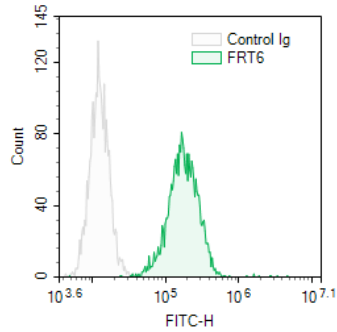

cell line:PC9  
2<sup>nd</sup> antibody:anti-mouse IgG FITC

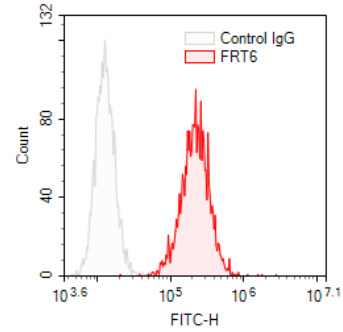

cell line:A375  
2<sup>nd</sup> antibody:anti-mouse IgG FITC

**Applications :** Flow Cytometry (FC), Immunoprecipitation (IP)

**Research Use :** For research use only, not for use in diagnostic procedures.

**Remark :** FC recommended concentration

(cell conc.  $5 \times 10^5$  cell / test & antibody conc.  $0.5 \mu\text{g}$  / test)

**Reference :**