

## *Anti-CD9 Monoclonal Antibody*

### *Purified antibody*

*(TSPAN-29)*

#### Product Information:

---

1) <b>Catalog #:</b>	MSFR106630
<b>Item Name:</b>	CD9 mAb FRT1(purified) 100ug
<b>Size :</b>	100µg (affinity-purified, protein G)
<b>Species :</b>	Mouse
<b>Product Code :</b>	CD9-Mo-MY001
<b>RRID :</b>	AB_2747365

---

**Clone name :** FRT1

**Other name :** MIC3; MRP-1; BTCC-1; DRAP-27; TSPAN29; TSPAN-29

**Formulation :** liquid. ( 0.93 mg / ml in PBS )

**Storage :** Store at 4°C. This antibody can also be aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thaw cycles. Non-hazardous. No MSDS required.


**Species :** mouse, monoclonal

**Isotype :** IgG1, κ


**Antigen :** MIA-PaCa2/Human

**Specificity :** human (others not tested)

This antibody is specific to CD9.

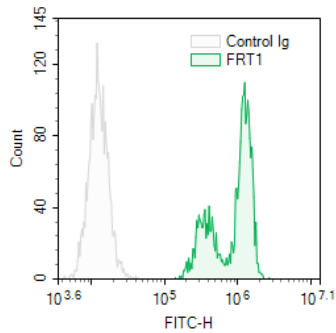

cell line:PC9  
2<sup>nd</sup> antibody:anti-mouse IgG FITC

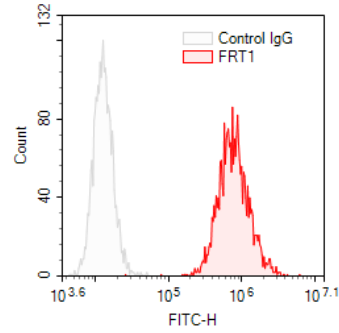

cell line:A375  
2<sup>nd</sup> antibody:anti-mouse IgG FITC

**Applications :** Flow Cytometry (FC), Immunoprecipitation (IP)

**Research Use :** For research use only, not for use in diagnostic procedures.

**Remark :** FC recommended concentration

(cell conc.  $5 \times 10^5$  cell / test & antibody conc.  $0.5 \mu\text{g}$  / test)

**Reference :**