

Anti-FAK

(focal adhesion kinase)

Code Number : FAK-Rb-Af340 (rabbit, RRID :)

FAK-GP-Af500 (guinea pig, RRID :)

Size : Size : 20 µg and 50 µg / See label on vial

(affinity-purified with antigen polypeptide)

Formulation : Liquid ; 200 µg/ml in PBS with 0.05% NaN₃.

Storage : Store at cool temp. (2-10°C)


The antibody can be stored at 2-10°C. The antibody can be also aliquotted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardrous. No MSDS required.

Species : rabbit / guinea pig, polyclonal

Antigen : rat, C-terminus 23 aa (AF020777)

Specificity : mouse (others not tested)

The specificity is shown by negative labeling in FAK-knockout brains.

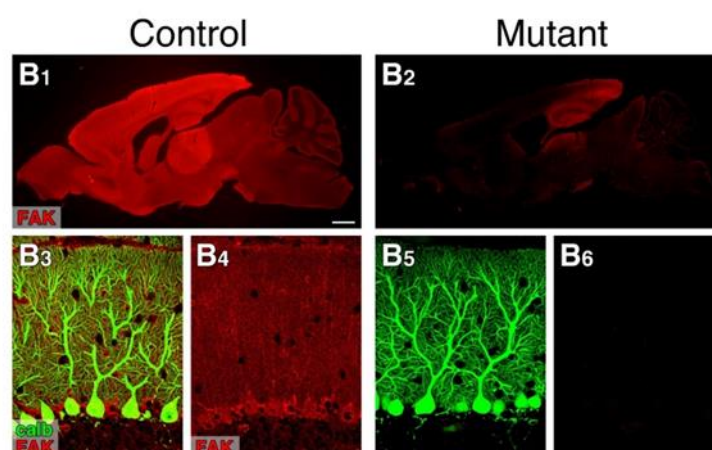

Applications : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : [NITTOBO MEDICAL CO., LTD](http://www.nittobogrp.com) Manufacture : FRONTIER INSTITUTE CO., LTD
TEL : +81-(0)3-4582-5451 email : nmd-tokuhan@nittobogrp.com

Research Use : For research use only, not for use in diagnostic procedures.

Remarks :

Reference : 1. Watanabe F, Miyazaki T, Takeuchi T, Fukaya M, Nomura T, Noguchi S, Mori H, Sakimura K, Watanabe M, Mishina M. (2008) Effects of FAK ablation on cerebellar foliation, Bergmann glia positioning and climbing fiber territory on Purkinje cells. **Eur J Neurosci.** 27:836-854.

ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : [NITTOBO MEDICAL CO., LTD](http://www.nittobo.com) Manufacture : FRONTIER INSTITUTE CO., LTD
TEL : +81-(0)3-4582-5451 email : nmd-tokuhan@nittobogrp.com