

Anti-FAK

(focal adhesion kinase)

Product Information:

1) Catalog #:	MSFR107020
Item Name:	FAK (focal adhesion kinase) pAb (Rb)
Size :	20µg (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	FAK-Rb-Af340
RRID :	AB_2858207
2) Catalog #:	MSFR107030
Item Name:	FAK (focal adhesion kinase) pAb (Rb)
Size :	50µg (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	FAK-Rb-Af340
RRID :	AB_2858207
3) Catalog #:	MSFR107040
Item Name:	FAK (focal adhesion kinase) pAb (GP)
Size :	20µg (affinity-purified with antigen polypeptide)
Species :	Guinea pig
Product Code :	FAK-GP-AF500
RRID :	AB_2858208
4) Catalog #:	MSFR107050
Item Name:	FAK (focal adhesion kinase) pAb (GP)
Size :	50µg (affinity-purified with antigen polypeptide)
Species :	Guinea pig
Product Code :	FAK-GP-AF500
RRID :	AB_2858208

Formulation: Liquid; 200 µg/ml in PBS with 0.05% NaN₃.

Storage: Store at cool temperature. (2-10°C)

The antibody can be stored at 2-10°C. The antibody can be also aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.


Species: rabbit / guinea pig, polyclonal

Antigen: rat, C-terminus 23 aa (AF020777)

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Specificity: mouse (others not tested)

The specificity is shown by negative labeling in FAK-knockout brains.

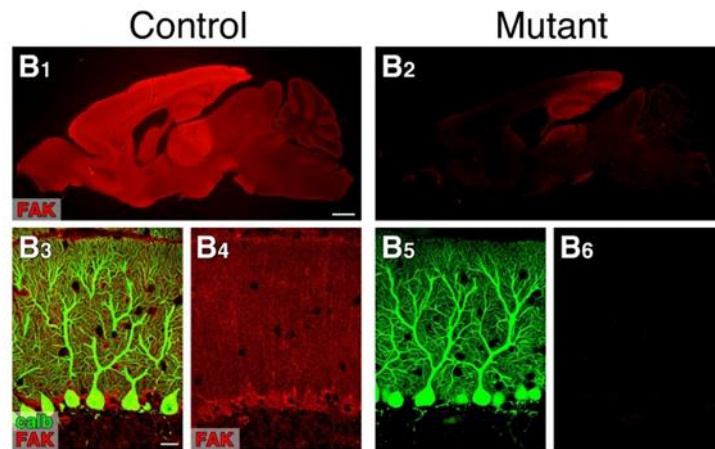

Applications: In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

Research Use: For research use only, not for use in diagnostic procedures.

Remarks:

Reference: 1. Watanabe F, Miyazaki T, Takeuchi T, Fukaya M, Nomura T, Noguchi S, Mori H, Sakimura K, Watanabe M, Mishina M. (2008) Effects of FAK ablation on cerebellar foliation, Bergmann glia positioning and climbing fiber territory on Purkinje cells. **Eur J Neurosci.** 27:836-854.

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD Manufacture : FRONTIER INSTITUTE CO., LTD
TEL : +81-(0)3-4582-5451 email : nmd-tokuhan@nittobogrp.com