

Anti-GST tag

Code Number : GSTAB

Size : 100 μ g (Affinity purified)


Formulation : PBS solution with 0.03% ProClin ® as preservative

Storage : Store at 4 °C. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.

Species : rabbit, polyclonal,

Antigen : GST tagged recombinant protein

Specificity : The antibody reacts with GST fusion protein

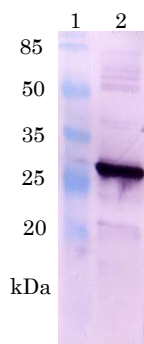

1, Molecular weight marker (kDa)

2, Affinity purified GST protein expressed in *E.coli*

Applications : Western blotting (x1/500 – x1/3000, Optimal dilution/concentration should be determined by the researcher.)

Not tested in other applications.

Research Use : For research use only, not for use in diagnostic procedures.

ご注意 : 本商品には 0.03 % の ProClin が入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD