

Anti-GluD2-N (GluR δ 2)

(glutamate receptor subunit δ 2)

Product Information:

1) Catalog #:	MSFR102630
Item Name:	GluD2-N (GluR δ 2) pAb (Rb) 20ug
Size :	20 μ g (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	GluRd2N-Rb-Af1460
RRID :	AB_2571760
2) Catalog #:	MSFR102640
Item Name:	GluD2-N (GluR δ 2) pAb (Rb) 50ug
Size :	50 μ g (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	GluRd2N-Rb-Af1460
RRID :	AB_2571760

Formulation : Liquid ; 200 μ g/ml in PBS with 0.05% NaN₃.

Storage : Store at cool temp. (2-10°C)

The antibody can be stored at 2-10°C. The antibody can be also aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.


Species : rabbit, polyclonal

Antigen : mouse GluD2,
N-terminal (extracellular) epitope,521-530aa (D13266)

Specificity : mouse (others not tested)

The specificity has been verified by knockout mouse brains.

Applications : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

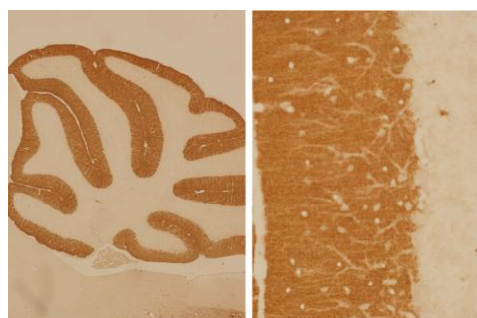

Research Use : For research use only, not for use in diagnostic procedures.

Remarks : For immunohistochemistry for neuronal iGluRs, users should adopt postembedding immunogold for electron microscopic detection and protease predigestion for light microscopic

ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。 誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

detection (see the below reference).

Reference :

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD