

## *Anti-GluD2-N (GluR $\delta$ 2)*

*(glutamate receptor subunit  $\delta$ 2)*

**Code Number** : GluRd2N-Rb-Af1460 (rabbit, RRID : AB\_2571760)

**Size** : 20  $\mu$ g and 50  $\mu$ g / See label on vial  
(affinity-purified with antigen polypeptide)


**Formulation** : Liquid ; 200  $\mu$ g/ml in PBS with 0.05% NaN<sub>3</sub>.

**Storage** : Store at cool temp. (2-10°C)

The antibody can be stored at 2-10°C. The antibody can be also aliquotted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardrous. No MSDS required.

**Species** : rabbit, polyclonal

**Antigen** : mouse GluD2,  
N-terminal (extracellular) eiptope,  
521-530aa (D13266)

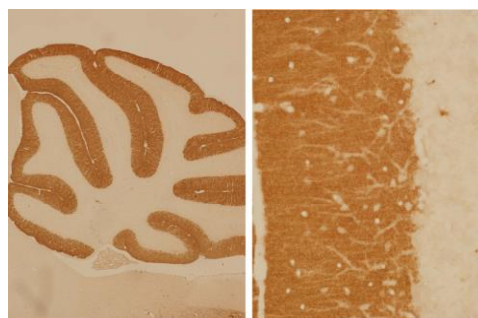

**Specificity** : mouse (others not tested)

The specificity has been verified by knockout mouse brains.

**Applications** : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

**Research Use** : For research use only, not for use in diagnostic procedures.

**Remarks** : For immunohistochemistry for neuronal iGluRs, users should adopt postembedding immunogold for electron microscopic detection and protease predigestion for light microscopic detection (see the below reference).

**Reference** :

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。