

Anti-GlyRa1

(glycine receptor $\alpha 1$ subunit)

Product Information:

| | |
|----------------|---|
| 1) Catalog #: | MSFR103030 |
| Item Name: | GlyRa1 pAb (Rb) 20ug |
| Size : | 20 μ g (affinity-purified with antigen polypeptide) |
| Species : | Rabbit |
| Product Code : | GlyRa1-Rb-Af840 |
| RRID : | AB_2571772 |
| 2) Catalog #: | MSFR103040 |
| Item Name: | GlyRa1 pAb (Rb) 50ug |
| Size : | 50 μ g (affinity-purified with antigen polypeptide) |
| Species : | Rabbit |
| Product Code : | GlyRa1-Rb-Af840 |
| RRID : | AB_2571772 |

Formulation : Liquid ; 200 μ g/ml in PBS with 0.05% NaN₃.

(affinity-purified with antigen polypeptide)

Storage : Store at cool temp. (2-10°C)



The antibody can be stored at 2-10°C. The antibody can be also aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.

Species : rabbit, polyclonal

Antigen : mouse GlyRa α 1, 380-403 aa (NM_020492)

Specificity : mouse (others not tested)

This is **selective to GlyRa α 1**, but not other GlyRa α subunits. This stains neuronal surface in the lower brainstem and spinal cord.

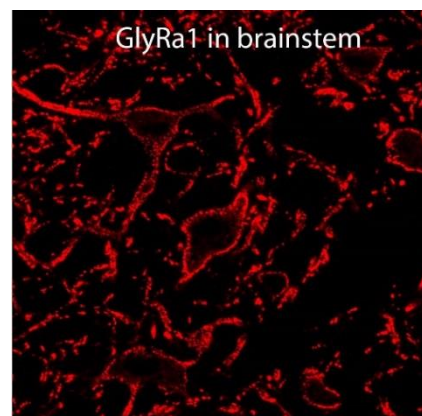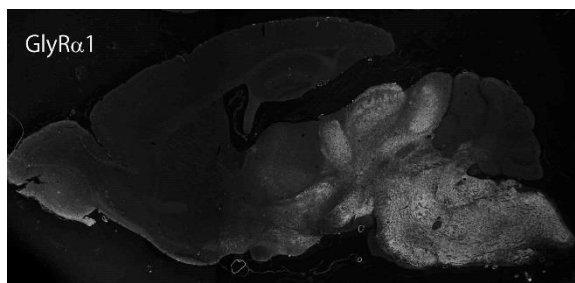

ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。 誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD

Applications : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

Research Use : For research use only, not for use in diagnostic procedures.

Remarks :

Reference :

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD