

Anti-Kv1.2

(Shaker-type potassium channel, type 2)

Product Information:

1) Catalog #:	MSFR103650
Item Name:	Kv1.2 pAb (Rb) 20ug
Size :	20µg (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	Kv1.2-Rb-Af240
RRID :	AB_2571789
2) Catalog #:	MSFR103660
Item Name:	Kv1.2 pAb (Rb) 50ug
Size :	50µg (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	Kv1.2-Rb-Af240
RRID :	AB_2571789

Formulation : Liquid ; 200 µg/ml in PBS with 0.05% NaN₃.

Storage : Store at cool temp. (2-10°C)

The antibody can be stored at 2-10°C. The antibody can be also aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.


Species : rabbit , polyclonal

Antigen : mouse, 482-496 aa (NM_008417)

Specificity : mouse (others not tested)

Immunoblot detects a protein band at 90 kDa In mouse, brains.

Applications : This can be used for immunohistochemistry and immunoblot.

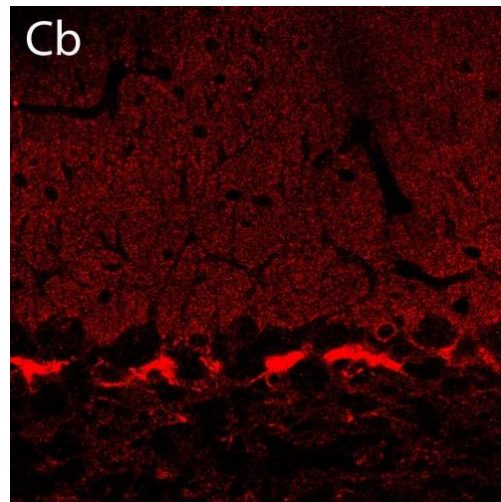

Research Use : For research use only, not for use in diagnostic procedures.

Remarks :

Reference : Iwakura A, Uchigashima M, Miyazaki T, Yamasaki M, Watanabe M: Lack of molecular-anatomical evidence for GABAergic influence upon axon initial segment of cerebellar Purkinje cells by the pinceau formation. **J. Neurosci.** 32:9438-9448.

ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。 誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD Manufacture : FRONTIER INSTITUTE CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com