

Anti-MT 3
(Metallothionein III)

Code Number : MT3-Mo-Af0218

Size : 20µg or 50µg (ProteinG purified)

Formulation : Liquid (50% Glycerol/PBS with 0.03% ProClin);

Storage : Store at 4 °C. The antibody can be stored at 4 °C. The antibody can be also aliquotted and stored at -80 °C for long-term storage. Avoid repeated freeze-thawing. Non-hazardrous. No MSDS required.


Species : mouse, monoclonal

Antigen :


recombinant human metallothionein 3

Specificity : Mouse/Human/Rat

This selectively stains liver, brain, and kidney.

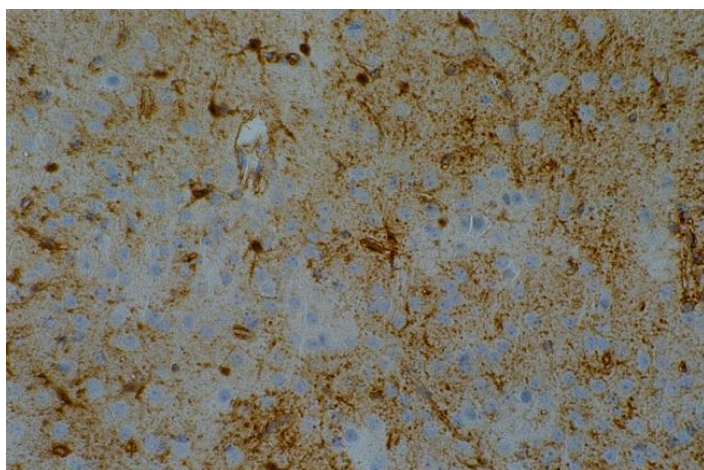

Mouse brain stained at 1µg/mL.

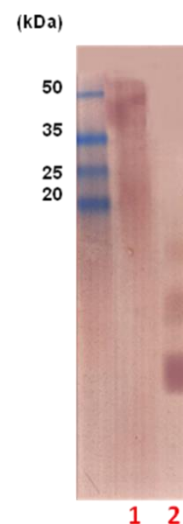

Lane1: MT-3 Apo, Lane2: MT-3 Zn

ご注意 : 本商品には 0.03 % の ProClin が入っています。 誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD

Applications: In general, protein purified antibody is used at 1 microgram/mL for immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

This antibody is used at 1.2microgram/mL for western blotting. Western blot analysis of MT-III. 2ug of MT-III proteins were electrophoresed on 15% gel NATIVE-PAGE, transferred to nitorocellulose membrane and then immunostained with this antibody.

Research Use: For research use only, not for use in diagnostic procedures.

Reference:

ご注意：本商品には 0.03 % の ProClin が入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD