

Anti-Metallothionein III

(MT3)

Product Information:

1) Catalog #:	MSFR077670
Item Name:	Metallothionein III mAb 100 μ g
Size :	100 μ g (Protein G purified)
Species :	Mouse
Product Code :	MT-3-MO-PG
RRID :	AB_2572280
2) Catalog #:	MSFR077680
Item Name:	Metallothionein III mAb 50 μ g
Size :	50 μ g (Protein G purified)
Species :	Mouse
Product Code :	MT-3-MO-PG
RRID :	AB_2572280

Formulation: Liquid (50% Glycerol/PBS with 0.03% ProClin);

Storage: Store at cool temperature. (2-10 °C)

The antibody can be stored at 2-10 °C. The antibody can be also aliquoted and stored at -80 °C for long-term storage. Avoid repeated freeze thawing. Non-hazardous. No MSDS required.

Species: mouse, monoclonal

Antigen:

recombinant human metallothionein 3


Specificity: Mouse/Human/Rat

This selectively stains liver, brain, and kidney.


ご注意：本商品には 0.03 % の ProClin が入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD

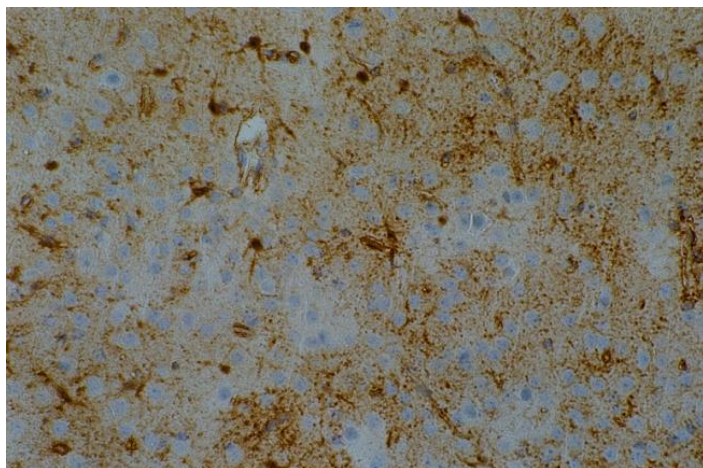

Mouse brain stained at 1µg/mL.

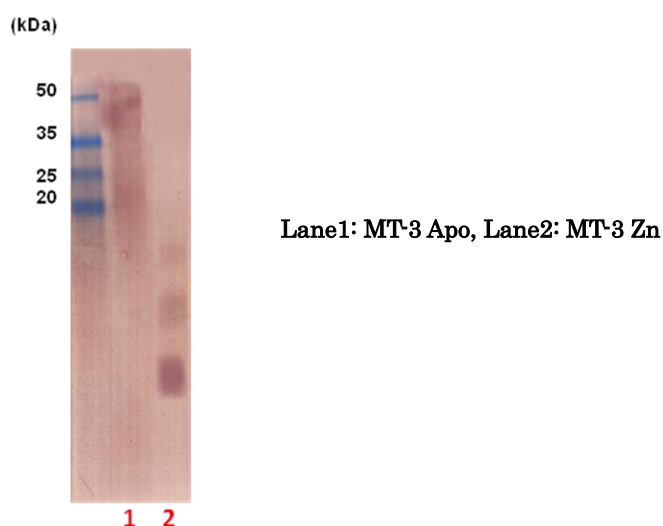

Applications: In general, protein purified antibody is used at 1 microgram/mL for immunohistochemistry. Users should determine the most appropriate concentration, because it depends on contents in given cells, tissues and organs.

This antibody is used at 1.2microgram/mL for western blotting. Western blot analysis of MT-III. 2ug of MT-III proteins were electrophoresed on 15% gel NATIVE-PAGE, transferred to nitrocellulose membrane and then immune-stained with this antibody.

Research Use: For research use only, not for use in diagnostic procedures.

Reference:

ご注意：本商品には 0.03 % の ProClin が入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD