

Anti-NPY (neuropeptide Y)

1) Catalog #:	MSFR104610
Item Name:	NPY (neuropeptide Y) pAb (Rb) 20ug
Size :	20µg (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	NPY-Rb-Af610
RRID :	AB_2571860
2) Catalog #:	MSFR104620
Item Name:	NPY (neuropeptide Y) pAb (Rb) 50ug
Size :	50µg (affinity-purified with antigen polypeptide)
Species :	Rabbit
Product Code :	NPY-Rb-Af610
RRID :	AB_2571860

Formulation : Liquid ; 200 µg/ml in PBS with 0.05% NaN₃.

Storage : Store at cool temp (2-10°C).

The antibody can be stored at 2-10°C. The antibody can be also aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.

Species : rabbit, polyclonal

Antigen : mouse NPY, 29-64 aa (NM_023456)


Specificity : This antibody stains cell bodies and fibers of NPY neurons in wild-type mice, but not in NPY-knockout mice.

Applications : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

Research Use : For research use only, not for use in diagnostic procedures.

Remarks :

Reference : 1) paper in submission

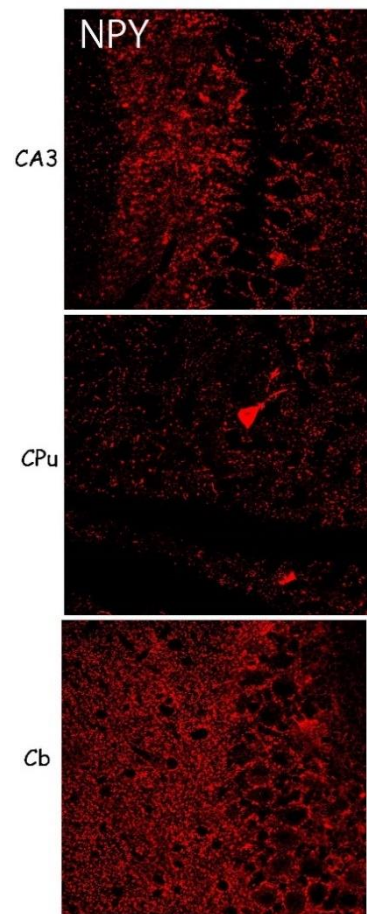

ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。