

Anti-POMC

(*pro-opiomelanocortin*)

Code Number : POMC-Rb-Af440 (rabbit, RRID :)

Size : 20 µg and 50 µg / See label on vial
(affinity-purified with antigen polypeptide)

Formulation : Liquid ; 200 µg/ml in PBS with 0.05% NaN₃.

Storage : Store at cool temp (2-10°C).


The antibody can be stored at 2-10°C. The antibody can be also aliquotted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardrous. No MSDS required.

Species : rabbit, polyclonal

Antigen : mouse, 27-76 aa (NM_001278581)

Specificity : mouse (others not tested)

This antibody stain cell bodies of the arcuate nucleus, and also label nerve fibers in the dorsomedial hypothalamic nucleus in the hypothalamus, paraventricular thalamic nucleus, and bed nucleus of the stria terminalis (intraamygdala division). The specificity is also supported by mutually exclusive labeling of cell bodies and fibers for POMC and agouti-related protein (AgRP) in these regions.

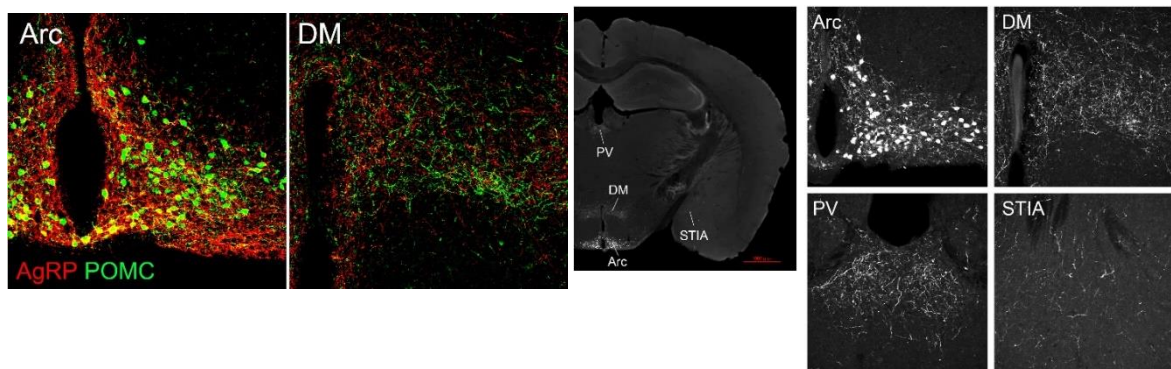

Applications : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : [NITTOBO MEDICAL CO., LTD](http://www.nittobomedical.com)
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD

Research Use : For research use only, not for use in diagnostic procedures.

Remarks :

Reference :

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : [NITTOBO MEDICAL CO., LTD](http://www.nittobomedical.com)
TEL : +81-(0)3-4582-5451
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD