

## *Anti-Tubulin $\beta$ -III*

### Product Information:

1)	<b>Catalog #:</b>	MSFR105960
	<b>Item Name:</b>	Tubulin $\beta$ III pAb (Rb) 20ug
	<b>Size :</b>	20 $\mu$ g (affinity-purified with antigen polypeptide)
	<b>Species :</b>	Rabbit
	<b>Product Code :</b>	TubIII-Rb-Af220
	<b>RRID :</b>	AB_2571848
2)	<b>Catalog #:</b>	MSFR105970
	<b>Item Name:</b>	Tubulin $\beta$ III pAb (Rb) 50ug
	<b>Size :</b>	50 $\mu$ g (affinity-purified with antigen polypeptide)
	<b>Species :</b>	Rabbit
	<b>Product Code :</b>	TubIII-Rb-Af220
	<b>RRID :</b>	AB_2571848
3)	<b>Catalog #:</b>	MSFR105980
	<b>Item Name:</b>	Tubulin $\beta$ III pAb (GP) 20ug
	<b>Size :</b>	20 $\mu$ g (affinity-purified with antigen polypeptide)
	<b>Species :</b>	Guinea pig
	<b>Product Code :</b>	TubIII-GP-Af760
	<b>RRID :</b>	AB_2571849
4)	<b>Catalog #:</b>	MSFR105990
	<b>Item Name:</b>	Tubulin $\beta$ III pAb (GP) 50ug
	<b>Size :</b>	50 $\mu$ g (affinity-purified with antigen polypeptide)
	<b>Species :</b>	Guinea pig
	<b>Product Code :</b>	TubIII-GP-Af760
	<b>RRID :</b>	AB_2571849

**Formulation** : Liquid ; 200  $\mu$ g/ml in PBS with 0.05% NaN<sub>3</sub>.

**Storage** : Store at cool temp (2-10°C).

The antibody can be stored at 2-10°C. The antibody can be also aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.


**Species** : rabbit / guinea pig, polyclonal

**Antigen** : mouse TubIII, C-terminal 17 aa  
(BC088749)

**Specificity** : mouse (others not tested)

This selectively stains neuronal dendrites, perikarya, and axons.

Tubulin beta-III C-17(rabbit)

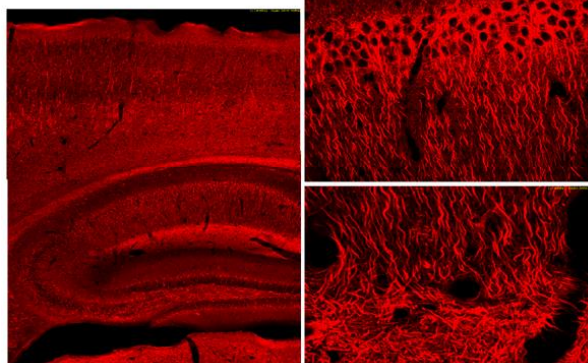

ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

**Applications** : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

**Research Use** : For research use only, not for use in diagnostic procedures.

**Remarks** : Rabbit and guinea pig antibodies are similar in the titer and specificity.

**Reference** :

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD  
TEL : +81-(0)3-4582-5451  
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD