

## *Anti-VR1*

(TRPV1, vanilloid receptor 1)

### Product Information:

---

1) <b>Catalog #:</b>	MSFR107180
<b>Item Name:</b>	VR1 (TRPV1) pAb (Rb)
<b>Size :</b>	20µg (affinity-purified with antigen polypeptide)
<b>Species :</b>	Rabbit
<b>Product Code :</b>	VR1-Rb-Af500
<b>RRID :</b>	AB_2858221
2) <b>Catalog #:</b>	MSFR107190
<b>Item Name:</b>	VR1 (TRPV1) pAb (Rb)
<b>Size :</b>	50µg (affinity-purified with antigen polypeptide)
<b>Species :</b>	Rabbit
<b>Product Code :</b>	VR1-Rb-Af500
<b>RRID :</b>	AB_2858221
3) <b>Catalog #:</b>	MSFR107200
<b>Item Name:</b>	VR1 (TRPV1) pAb (GP)
<b>Size :</b>	20µg (affinity-purified with antigen polypeptide)
<b>Species :</b>	Guinea pig
<b>Product Code :</b>	VR1-GP-Af590
<b>RRID :</b>	AB_2858222
4) <b>Catalog #:</b>	MSFR107210
<b>Item Name:</b>	VR1 (TRPV1) pAb (GP)
<b>Size :</b>	50µg (affinity-purified with antigen polypeptide)
<b>Species :</b>	Guinea pig
<b>Product Code :</b>	VR1-GP-Af590
<b>RRID :</b>	AB_2858222

---

**Formulation** : Liquid ; 200 µg/ml in PBS with 0.05% NaN<sub>3</sub>.

**Storage** : Store at cool temp (2-10°C).

The antibody can be stored at 2-10°C. The antibody can be also aliquoted and stored at -80°C for long-term storage. Avoid repeated freeze-thawing. Non-hazardous. No MSDS required.

**Species** : rabbit / guinea pig, polyclonal

**Antigen** : mouse, C-terminal 15 aa (EDAEVFKDSMAPGEK, NM001001445)

**Specificity** : mouse (others not tested)

This antibody strongly stains cell bodies and axons of primary sensory neurons in the dorsal root ganglion and spinal dorsal horn, respectively.



ご注意 : 本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

**Applications** : In general, affinity-purified antibody is used at around 1 microgram/ml for immunoblot and immunohistochemistry. The most appropriate concentration should be determined by users, because it depends on contents in given cells, tissues and organs.

**Research Use** : For research use only, not for use in diagnostic procedures.

**Remarks** :

**Reference** :

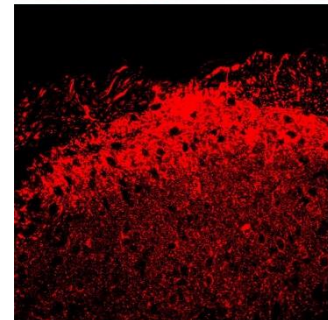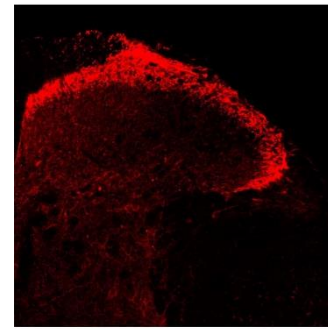

ご注意：本商品には 0.1 % 未満のアジ化ナトリウムが入っています。誤って目や口に入ったり、皮膚に付着した場合は大量の水で洗い流してください。

Distributor : NITTOBO MEDICAL CO., LTD  
TEL : +81-(0)3-4582-5451  
email : nmd-tokuhan@nittobogrp.com

Manufacture : FRONTIER INSTITUTE CO., LTD